CE₀₁₂₃

Metacross™RX

PTA Balloon Dilatation Catheter
Cathéter de dilatation à ballonnet pour PTA
PTA-Balloon dilationskatheter
Catéter de dilatación de balón para ATP
Cateter de Dilatação de Balão para ATP
Catetere con palloncino di dilatazione per PTA
PTA ballondilatatiekatheter
PTA ballongdilations kateter
PTA ballon-dilationskateter
PTA ballongdilatasjonskateter
PTA-pallolaajennuskatetri
Καθετήρας Διαστολής PTA με Μπαλόνι
Периферический баллонный катетер ЧТА
Balonowy cewnik dylacyjny do PTA
PTA Tágító Ballonkatéter
Balónkový dilatační katétr pro PTA
PTA balónikový dilatačný katéter
PTA Balon Dilatasyon Kateteri
PTA balloondilatatsiooni kateeter
PTA dilatācijas balonkatetrs
PTA kateteris su išplečiamuoju balionėliu
Balonski kateter za dilataciju PTA
PTA balon dilatacion kateter
Cateter de dilatare pe balon pentru PTA
PTA дилатационен балон-катетър
Дилатаційний балон-катетер для ЧТА
Balonski dilatacijski kateter za PTA

**TERUMO**

September, 2021

PM-05580
1/115 (3209-3)

REF	SN	LOT	
Catalogue number Référence du catalogue Bestellnummer Número de catálogo Referência Codice prodotto Catalogusnummer Artikelnnummer Varenummer Artikkelnummer Tuotekoodi Κωδικός Προϊόντος Каталожный номер Kod produktu Termékkód Kód výrobku Kód výrobku Ürün kodu Kood Koda Nr. Katalogo Nr. Koda izdelka Šifra proizvoda Codul produsului Каталожен номер Номер за каталогом Kataloški broj	Serial number Numéro de série Seriennummer Número de serie Número de série Numero di serie Serienummer Serienummer Serienummer Serienummer Sarjanumero Αριθμός παραγωγής Серийный номер Numer serii Sorozatszám Sériové číslo Sériové číslo výrobku Seri Numarası Seerianumber Sērijas Nr. Serijos Nr. Serijska številka Serijski broj Număr de serie. Сериен номер Серийный номер Serijski broj	Lot number Numéro de lot Chargenbezeichnung Número de lote Número de lote Numero di lotto Lotnummer Batchkod Batchnummer Batch number Eränumero Αριθμός παρτίδας Номер партии Numer serii Tételszám Číslo výrobní serie Identifikačné číslo Lot numarası Partii nr. Sērijas Nr. Partijos kodas Številka serije Serijska oznaka Număr lotului Партиден номер Код партії Šifra serije	Use by A utiliser jusque Verwendbar bis Caducidad Data limite de uso Utilizzare entro Gebruiken voor Användes före Udløbsdato Brukes før-dato Eräpäivä Ημερομηνία Λήξης Годен до Użyć przed datą Felhasználható Použit do Použit do Son kullanim tarihi Säilivusaeg Derīguma termiņš Sunaudoti iki Uporabno do Upotrebljivo do Data de expirare Срок на годност Використати до Upotrijebiti do
STERILEEO			P
Sterilized using ethylene oxide Stérilisé à l'oxyde d'éthylène Sterilisation mit Ethylenoxid Esterilizado con óxido de etileno Esterilizado com oxido de etileno Sterilizato con ossido di etilene Gesteriliseerd met ethyleenoxide Steriliserad med etylenoxid Ætylenoxid steriliseret Sterilisert med etylen oksid Steriloitu etyleenoksidilla Αποστειρωμένο με αιθυλαινοξειδιο Стерилизовано окисью этилена Sterylizowane tlenkiem etylenu Etilénoxidál sterilizálva Sterilizováno ethylenoxidem Sterilizované etylénoxidom Etilen Oksitle sterilize edilmiştir Steriliseeritud etüleenoksiidiga Sterilizēts ar etilēna oksīdu Sterilizuota naudojant etileno oksidą Sterilizirano z etilen oksidom Sterilisano etilen oksidom Sterilizat cu oxid de etilena Стерилизовано с этилен оксид Простерилизовано за допомогою оксиду етилену Sterilizirano etilen-oksidadom		Do not reuse Ne pas réutiliser Nicht zur Wiederverwendung No reutilizar Não reutilizar Monouso Niet opnieuw gebruiken Ingen återanvändning Må ikke genbruges Må ikke gjenbrukes Kertakäyttöinen Μόνο για μια χρήση Не использовать повторно Nie używać powtórnie Csak egyszeri használatra Nelze použiť opakovaně Na jednorázové použitie Tek kullanımlık Ainult ühekordseks kasutuseks Tikai vienreizējai lietošanai Pakartotina nenaudoti Samo za enkratno uporabo Samo za jednokratnu upotrebu Doar pentru unica folosinta Само за еднократна употреба Використовувати повторно заборонено Ne upotrebljavati ponovno	Pressure Pression Druck Presión Pressão Pressione Druk Tr yck Tr yk Trykk Paine Πίεση Давление Ciśnienie Nyomás Tlak Tlak Basiņš Röhk Spiediens Slegis Tlak Pritisak Presiune Налягане Тиск Tlak

	
Consult instructions for use Lire le mode d'emploi Gebrauchsanweisung beachten Consulte las instrucciones de uso Consultar instruções de uso Consultare le istruzioni per l'uso De gebruiksaanwijzing raadplegen Läs bruksanvisningen Læs brugsanvisningen Les bruksanvisning før bruk Lue käyttöohjeet Διαβάστε τις οδηγίες χρήσης Прочитайте инструкции по применению Sprawdź w instrukcji stosowania Olvassa el a használati útmutatót! Postupujte podle pokynů Prečítajte si návod na použitie Kullanma talimatına başvurunuz Tutvu kasutusjuhendiga. Pirms lietošanas izlasiet instrukciju Skaitykite naudojimo instrukciją Glej navodila za uporabo Pogledaj instrukcije za upotrebu A se consulta instrucțiunile de folosință Прочетете указанията за употреба Ознайомтеся з інструкціями для застосування Pogledati upute za upotrebu	Manufacturer Fabricant Hersteller Fabricante Fabricante Fabbicante Fabrikant Tilverkare Produceret af Produsent Valmistaja Κατασκευαστής Изготовлено Wytwórca Gyártó Výrobeno Výrobca Üretici firma Toodetud Razõts Gamintojas Izdelevalac Proizvodjač Fabricat de Произведено от Виробник Proizvođač
NP	RBP
Nominal Pressure Pression nominale Nenndruck Presión nominal Pressão nominal Pressione nominale Nominale druk Nominellt tryck Nominelt tryk Nominelt trykk Nominaalipaine Ονομαστική Πίεση Номинальное давление Ciśnienie znamionowe Névleges nyomás Jmenovitý tlak Menovitý tlak Nominal Basiņ Nominaalrõhk Nominālais spiediens Minimalus slėgis Nominalni tlak Nominalni pritisak Presiune nominală Номинално налягане Номинальний тиск Nominalni tlak	Rated Burst Pressure Pression nominale de rupture Zerreißdruck Presión de rotura Pressão nominal de ruptura Pressione massima di sicurezza Bepaalde scheurdruk Nominellt bristningstryck Nominel øvre trykgrænse Nominelt sprengtrykk Nimellinen rikkoutumispaine Ονομαστική Πίεση Ρήξης Расчетное давление разрыва Znamionowe ciśnienie rozrywające Becsült szakítónyomás Jmenovitý tlak prasknutí Menovitý tlak prasknutia Anma Patlama Basıncı Arvestuslik purunemisrõhk Aprēķinātais sprādziena spiediens Minimalus išplėtimo slėgis Naznačen delovni tlak Granični pritisak pucanja Presiunea nominală de explozie Номинално налягане на разкъсване Номинальний тиск розриву Nazivni tlak prsnuća

 Contents	
Contents Contenu Stück Unidades Conteúdo Unità Stuks Styck Indhold Innhold Sisältö Περιεχόμενα Содержание Zawartość Tartalom Obsah Obsah İçindekiler Sisu Satus Sudėtis Vsebina Sadržaj Conținut Зміст Вміст Sadržaj	Do not use if package is damaged Ne pas utiliser si l'emballage est endommagé Inhalt beschädigter Packung nicht verwenden No usar si el paquete está dañado Não usar se a embalagem estiver danificada Non utilizzare se la confezione risulta danneggiata Niet gebruiken wanneer de verpakking beschadigd is Använd inte om förpackningen är skadad Må ikke anvendes, hvis pakningen er beskadiget Må ikke brukes hvis forpakningen er ødelagt Älä käyttää jos pakkaus on vahingoittunut Μην το χρησιμοποιήσετε εάν η συσκευασία είναι κατεστραμμένη Не использовать, если упаковка повреждена Nie stosować w przypadku uszkodzenia opakowania Ne használja, ha a csomagolás sérült! Nepoužívejte, bylo-li balení poškozeno Nepoužívejte, ak je obal poškodený Paket hasarlıysa kullanmayınız Årge kasutage, kui pakend on kahjustatud Nelietot, ja sterila isepakojums ir bojāts Jei pakuotė pažeista, nenaudoti Ne uporabljajte ga, če je ovojnina poškodovana. Ne upotrebljavati ukoliko je pakovanje oštećeno A nu se utiliza în cazul în care ambalajul este deteriorat Не използвайте, ако опаковка е повредена Не використовувати, якщо пакування пошкоджено Ne upotrebljavati ako je pakiranje oštećeno
	
Balloon diameter Diamètre du ballon Durchmesser des Ballons Diámetro del balón Diâmetro do balão Diametro del palloncino Ballondiameter Ballongdiameter Ballondiameter Ballongdiameter Pallon halkaisija Διάμετρος μπαλονιού Диаметр баллона Średnica balonu Ballon átmérő Průměr balónku Priemer balóna Balon çapı Ballooni diameeter Balona diametrs Ballono diametras Premier balona Prečnik balona Diametrul balonului Диаметър на балона Диаметр балону Promjer balona	Balloon length Longueur du ballon Länge des Ballons Longitud de balón Comprimento do Balão Lunghezza del palloncino Ballonlengte Ballonglängd Ballonlængde Ballonglengde Pallon pituus Μήκος Μπαλονιού Длина баллона Długość balonu Ballon hossz Délka Balonku Dĺžka balóna Balon Uzunluğu Ballooni pikkus Balona garums Ballono ilgis Dolžina balona Dužina balona Lungimea balonului дължина на балона довжина балону Duljina balona

EC REP	2 STERILIZE
Authorized Representative in the European Community Rep. CE (Mandataire européen) Bevollmächtigter in der Europäischen Gemeinschaft Representante autorizado de la Comunidad europea Representante autorizado na CE Rappresentante autorizzato nella comunità europea Gemachtigde vertegenwoordiger in de EU (Europese Unie) Auktoriserad återförsäljare inom EU-området Autoriseret forhandler indenfor EU området Autoriseret representant i Det europeiske Fællesmarked EY edustaja Εξουσιοδοτημένος αντιπρόσωπος στην Ευρωπαϊκή Κοινότητα Авторизованный представитель Европейского сообщества Autoryzowany przedstawiciel w Unii Europejskiej Az Európai Közösségben a képviseletre jogosult Zástupce Evropského společenství Predstaviteľ Európskej únie Yetkili AB temsilcisi Esindaja EL-s ES pārstāvis Įgaliotas atstovas Europos bendrijoje Zastopnik v Evropski Skupnosti Predstavnik u Evropskoj Uniji Reprezentant UE Упълномощен представител в ЕС Уповноважений представник в Европейському Союзу Ovlašteni predstavnik u europskoj zajednici	Do not resterilize Ne pas restériliser Nicht wieder sterilisieren No reesterilizar Não re-esterilizar Non risterilizzare Niet hersteriliseren Får ej omsteriliseras Må ikke reesteriliseres Må ikke re-steriliseres Ei uudelleen sterilointia Να μην επαναποστείρωθεί Не стерилизовать повторно Nie sterylizować ponownie Nem sterilizálható újra! Nesterilizujte opakovane Nesterilizujte Yeniden sterilize etmeyin Ärge steriliseerige uuesti Nesterilizēt atkārtoti Pakartotinai nesterilizuoti Ne reesterilizirati Ne resterilisati Nu resterilizați Да не се стерилизира повторно Повторно не стерилизувати Ne ponovo sterilizirati
	
Manufacturing Site Site de fabrication Herstellungsort Lugar de fabricación Site do fabricante Sito di produzione Productielocatie Tillverkningsplats Produktionssted Produksjonssted Valmistuspaikka Τόπος Κατασκευής Произведено Miejsce produkcji Gyártóhely Misto výroby Miesto výroby Üretim Yeri Tootmiskoht Ražotne Gambybos vieta Proizvodno mesto Mesto proizvodnje Locul de Productie Место на производство Виробнича дільниця Mjesto proizvodnje	Keep dry Conserver au sec Vor Nässe schützen Mantener seco Manter seco Conservare in luogo asciutto Droog bewaren Håll torrt Holdes tør Holdes tørt Pidä kuivassa Διατηρήστε το στεγνό Беречь от влаги Przechowywać w miejscu suchym Nedvességtől védett helyen tartandó Udržujte v suchu Udržujte ja suchom mieste Kuru tutunuz Hoida kuivas Sargāt no mitruma Laikyti sausoje vietoje Izogibajte se dežju Čuvati na suvom mestu A se păstra în locuri ferite de umezeală Съхранявайте на сухо Зберігати в сухому місці Držati suho

Keep away from sunlight
Craint la chaleur
Vor Sonnenlicht schützen
Mantener alejado de la luz solar
Manter ao resguardo da luz solar
Tenere lontano da fonti di calore
Niet blootstellen aan zonlicht
Aktas för solljus
Må ikke udsættes for sollys
Unngå sollys
Suojattava auringonvalolta
Φυλάξτε το μακριά από το ηλιακό φως
Избегать прямого солнечного света
Chronić przed bezpośrednim działaniem promieni słonecznych
Napfénytől védett helyen tartandó
Chraňte před slunečním zářením
Udržujte mimo slnečného žiarenia
Güneş ışığından uzak tutunuz
Kaitsta otsese päikesevalgus eest
Sargāt no saules stariem
Saugoti nuo saulės šviesos
Izogibajte se sončni svetlobi
Ne izlagati sunčevoy svetlosti
A se pästra la umbrä
Пазете от сплънчева светлина
Зберігати подаль від сонячного проміння
Drżati dalje od sunčevog svjetla

Rapid Exchange Type
Echange rapide
Rapid Exchange Typ
del Tipo de Intercambio Rápido
Tipo de Troca Rápida
Tipo a Scambio Rapido
Rapid Exchange Type
Snabbutbytestyp
Rapid Exchange Type
Hurtig Utskifting
Nopeasti Vaihdehtava
Ταχείας Ανταλλαγής
Быстрой Замены
Typ szybkiej wymiany
Gyorscserele rendszeru
Rychlou Výmenou
Rýchlo Vymeniteľný
Rapid Exchange Tipi
Kiirvahetustüüpi
Atras Apmainas
Greitai Keiciamas
Hitro Izmenjavo
Rapid Exchange Type
Schimb rapida
Бързосменния
Швидкої Заміни
Tip za brzu razmjenu

Rx ONLY

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
[PRODUCT DESCRIPTION]

Metacross RX is a Percutaneous Transluminal Angioplasty (PTA) Balloon Dilatation Catheter for peripheral indications. Metacross RX is a rapid exchange (RX) type and the maximum diameter of the compatible guidewire is 0.035 inches (0.89 mm).

The distal section of this catheter consists of a balloon and dual lumen shafts (co-axial). The outer lumen is used for inflation and deflation of the balloon with contrast medium diluted with saline solution. The inner lumen (the guidewire lumen; from the distal tip to the guidewire port) is for inserting a compatible guidewire to facilitate advancing the catheter through the stenotic lesion or stent to be dilated. The proximal section is a single lumen shaft with a single luer port hub for connecting an inflation/deflation device. Inside the balloon, two radiopaque markers, which indicate the working length of the balloon, are placed to guide the physician for positioning the balloon properly in the targeted lesion under fluoroscopy. Metacross RX is available in various balloon sizes by its diameter and length at recommended inflation pressures (Nominal pressure). Upon inflation, the balloon diameter varies according to the inflating pressure. The balloon compliance chart of Metacross RX is provided in Table 1 at the end of this Instructions for Use. In the sterilization package, the device is equipped with the protective materials, a balloon protective tube and a stylet, which are to be removed before use. The balloon protective tube is placed over the balloon folded in a low profile and the stylet is inserted into the guidewire lumen to prevent collapsing the catheter shaft. The catheter's distal tip is tapered to facilitate the catheter's advancing into the stenotic lesion or stent.

The flushing needle is provided as an accessory of this device for flushing and filling the guidewire lumen with heparinized saline.

<NAME OF EACH PART>

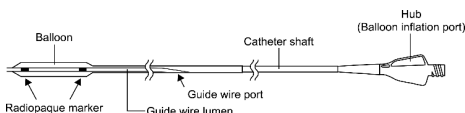Sterile and non pyrogenic in an undamaged and unopened blister pouch. This device is sterilized by ethylene oxide.

Recommended guidewire diameter:

Maximum diameter: 0.035"(0.89 mm)

Recommended inflation pressure and maximum inflation pressure

Nominal pressure:

- Balloon diameter: smaller than 6 mm: 10atm
- Balloon diameter: 6 mm and larger: 8atm

Rated burst pressure:

- Balloon diameter: smaller than 5 mm: 20atm
- Balloon diameter: 5 mm : 18atm
- Balloon diameter: 6 mm, 7 mm, 8 mm : 14atm
- Balloon diameter: 9 mm and larger: 12atm

[INDICATIONS]

The Metacross RX PTA balloon dilatation catheter is intended to dilate stenosis in the iliac, femoral, ilio-femoral, popliteal, infra-popliteal and renal arteries, and for the treatment of obstructive lesions of native or synthetic arteriovenous dialysis fistulae. This device is also indicated for stent post-dilatation in the peripheral vasculature.

Note: The verification test for stent post-dilatation of Metacross RX was conducted using Complete® SE Vascular Stent System (Medtronic, Inc), S.M.A.R.T.® CONTROL® Vascular Stent System (Cordis Corporation) and Epic™ Vascular Self-Expanding Stent System (Boston Scientific Corporation).

[CONTRA-INDICATIONS]

1. Patients who have developed anastomotic stenosis within one month after AV shunt construction.
2. Patients having lesions which communicate with a pseudoaneurysm.
3. Patients who are pregnant.
4. Patients who cannot tolerate antiplatelet therapy or anticoagulant therapy.

[WARNINGS AND PRECAUTIONS]

[Warnings]

1. This device is for single use only. Do not reuse. Do not resterilize. Do not reprocess. Reprocessing may compromise the sterility, the biocompatibility and the physical integrity of the device.
2. Do not use if the product or the unit packaging has been damaged.
3. Do not inflate the balloon exceeding the diameter of the blood vessel proximal and distal to the stenotic lesion.
4. Do not inflate the balloon to a pressure exceeding the rated burst pressure.
5. Do not use this catheter in the central circulatory system.

[Precautions related to procedures]

1. Do it carefully and slowly to remove this device out of the carrier tube and to remove protection parts from this device, Careless handling may damage the balloon and/or the catheter shaft of this device to impair balloon inflation and/or deflation.
2. Inflate the balloon with the contrast medium diluted with saline in the ratio of 1:1 (hereafter, inflation fluid). No gaseous medium such as air should be used for inflation. (Balloon may be inflated irregularly.)
3. Prior to use, expel all the air in the balloon and balloon inflation lumen and replace it with the inflation fluid. (In case of incomplete air removal, balloon inflation state can not be observed under fluoroscopy.)
4. In case of serious stenotic lesions such as calcified lesions, the blood vessel can not be dilated fully. Do not apply a pressure exceeding the rated burst pressure in such a case. (The balloon may burst and the debris may remain inside the body.)
5. As for the medical devices used in conjunction with this catheter, follow the instructions for use of such devices. (When the balloon catheter is used for post-dilatation of stent during peripheral stent placement procedure, refer to the manufacturer's instructions for use.)
6. For insertion into stent, withdrawal from stent or in-stent dilatation of the balloon, proceed with caution under fluoroscopy. (Doing so without fluoroscopy may damage this catheter or injure the blood vessel.)
7. Do not attempt to pass the Metacross RX PTA balloon dilatation catheter through a smaller sized introducer sheath than indicated on the label. Refer to product labeling.

[Precautions during usage]

1. Use the device immediately after the sterile package is opened.
2. The entire procedure should be carried out aseptically.
3. This catheter should be used only by physicians skilled in percutaneous vascular therapy.
4. Select the catheter with appropriate balloon size (balloon diameter, balloon length) according to the following (procedures) criteria.
 - 1) The inflated balloon diameter shall not exceed the inner vessel diameter proximal and distal to the lesion.
 - 2) The length of the inflated balloon shall not exceed the length of the lesion to be treated.

(Inappropriate catheter selection may lead to acute vascular occlusion.)

5. Use this catheter only when emergency surgery can be performed at any time.
6. The physician in charge of the procedure should determine the duration and number of balloon inflations based on his/her past experiences.
7. Heparinized saline should be infused for anti-coagulation while this balloon catheter is inserted in the blood vessel.
8. Never inflate the balloon or advance the guidewire with the catheter bent or kinked.
9. Manipulate the catheter carefully in the blood vessel verifying the location and movement of its tip under fluoroscopy or DSA (Digital Subtraction Angiography) monitoring.
10. Always inflate the balloon while observing its inflation state under fluoroscopy. Inflate the balloon carefully while verifying its movement with the radiopaque markers. If any abnormality is found such as a little resistance, immovable catheter during manipulation or kinked distal tip of the guidewire, the procedure should be discontinued immediately and this catheter should be removed slowly together with the guidewire.
11. Do not insert or remove the catheter rapidly. (Operating rapidly may damage the catheter or injure the vascular intima.)
12. Do not move the catheter with the balloon inflated. (Moving the catheter with the balloon inflated may result in balloon burst or catheter shaft breakage.)
13. Do not inflate or deflate the balloon rapidly in the blood vessel. (Rapid inflation or deflation may damage the blood vessel or cause the balloon to burst resulting in the debris left inside the body.)
14. Do not use an injector (automatic injection device) for inflating the balloon.
15. Precautions should be taken to prevent any damage to the catheter by a surgical knife or scissors.
16. If any abnormality such as strong resistance is experienced while manipulating the catheter, the procedure should be discontinued immediately. The cause should be verified and appropriate measures should be taken. (Continuing the operation with excessive force may result in damage to the catheter or in vascular wall injury.)
17. Always pay attention to the patient's condition and monitor the body temperature, pulse, and breathing. When any abnormality is found, discontinue the use of the catheter immediately or take appropriate measures for the patient's condition on the discretion of the physician.
18. After use dispose the catheter as medical waste according to hospital procedures.
19. Do not use catheters containing organic solvents or oleaginous contrast media. Contact with these agents may lead to damage of the catheter.
20. While manipulating this catheter, do not twist or rotate the catheter.
21. While inserting this device into the blood vessel or removing this device from the blood vessel, make sure that the balloon is completely deflated. (A device with larger and longer balloon requires a longer deflation time.)
22. If resistance is felt during post procedural withdrawal of this device, it is recommended to withdraw the entire system together with the introducer / guide sheath.
23. Always use the catheter along a guide wire inserted into the guide wire lumen.

[Operational Instructions]

Materials typically required for PTA with the Metacross RX PTA balloon dilatation catheter include:

- Guidewire(s) of appropriate diameter and length
- Appropriate introducer/guide sheath and dilator set
- Vial of contrast medium
- Vial of sterile saline
- Inflation device with manometer
- Luer-lock syringe
- Three-way stopcock

[Operation method or instructions for use]

1. Selection of the catheter
Prior to use, verify the blood vessel inner diameter proximal and distal to the stenotic lesion under fluoroscopy, and select the catheter with appropriate balloon size. If two sizes are applicable to the diameter verified, select the smaller size.
2. Preparations
 - 1) After aseptically removing the catheter from the package container, remove the protective materials from the catheter.
 - 2) Replace the air in the balloon and balloon inflation lumen with inflation fluid according to the following procedures.
 - a) Attach the inflation device (not included in this kit) filled with the inflation fluid onto the balloon inflation port.
 - b) After applying negative pressure for approximately 10 seconds to the inflation device, direct the catheter shaft tip downward and release the negative pressure slowly to allow the fluid to fill the balloon and balloon inflation lumen and to expel the air.
 - c) When any air is found in the balloon and balloon inflation lumen, repeat the procedure b) until the air is expelled completely.
 - d) Detach the inflation device from the balloon inflation port and remove air from the cylinder of the inflation device.
 - e) Reattach the inflation device onto the balloon inflation port and maintain negative pressure. Verify that the air no longer returns into the syringe.
 - 3) Flush the guidewire lumen with heparinized saline using the flushing needle and fill the lumen with the heparinized saline.
 - 4) Immerse the entire Metacross RX PTA balloon dilatation catheter in saline.
3. Insertion and inflation of the balloon catheter
Prior to use, please check if this catheter is compatible to the guide catheter (not included in this kit), the introducer sheath (not included in this kit), and the guidewire (not included in this kit) by the following table.

Catalogue number	Balloon diameter (mm)	Compatible sheath (Fr/mm)	Maximum diameter for the guidewire (inch/mm)
BD-L30020MR	3 mm, 4 mm, 5 mm, 6 mm, 7 mm, 8 mm	6 Fr / 2.00 mm	0.035 inch (0.89 mm)
BD-L70200ER	9 mm, 10 mm, 12 mm	7 Fr / 2.33 mm	

When using the catheter inserted in the introducer sheath or guide catheter

- 1) When using the catheter in the guide catheter, loosen the hemostatic valve while keeping it tight enough to prevent blood leakage.
- 2) Gradually insert the catheter, with the balloon deflated completely, into the introducer sheath or guide catheter over the guidewire and advance to the distal periphery of the lesion. Under fluoroscopy, advance the catheter slowly to make its radiopaque marker in the balloon reach the target site.

- 3) Under fluoroscopy, verify the position of the radiopaque marker in the balloon and determine the dilatation site.
- 4) Attach the inflation device securely to the balloon inflation port.
- 5) Inflate the balloon gradually to the target diameter verifying the inflation pressure and dilate the stenotic lesion (for inflated balloon diameters, consult the balloon compliance chart provided in Table 1. at the end of this Instructions for Use).
4. Removal of the catheter
 - 1) After applying negative pressure to deflate the balloon completely, remove the catheter.
5. Re-insertion of the catheter
 - 1) If a re-insertion of the removed catheter with deflated balloon is required, ensure the balloon is completely deflated and folded before re-inserting it into the patient's blood vessel. On re-insertion use the catheter as per "Insertion and inflation of the balloon catheter" section herein.
 - 2) If any resistance is encountered while re-inserting the catheter, stop advancing it and carefully remove it out of the patient's blood vessel. Prepare to use a new Metacross RX in accordance with the procedures provided above in this section of this Instruction for Use.

[Complications]

Device Failures:

- Balloon rupture
- Insufficient inflation / deflation of the balloon
- Breakage of the balloon and / or the catheter shaft
- Difficulty in removing the device
- Leakage of inflation fluid

Adverse Events:

- Local or systemic infections
- Local internal bleeding or hematoma
- Intimal rupture
- Vascular dissection
- Vascular perforation
- Vascular rupture
- Aneurysm
- Arrhythmia
- Acute vascular occlusion
- Venous thromboembolism
- Vasospasm
- Formation of pseudoaneurysm
- Arteriovenous fistula
- Bleeding requiring transfusion
- Allergic response to the contrast media / renal failure
- Ache or pressing pain
- Death
- Embolism
- Endocarditis
- Fever
- Hypertension / Hypotension
- Inflammation
- Myocardial Infarction
- Sepsis
- Shock
- Stroke
- Transient Ischemic Attack

[Storage and expiration date]

1. Store in a cool, dry, dark place. Avoid exposure to water, direct sunlight.
2. The expiration date is indicated on the box. Do not use after the expiration date.

[Package]

1 set/box

[Name and address of manufacturer]

Manufacturer : KANEKA CORPORATION

Address: 3-18, 2-chome, Nakanoshima, Kita-ku, Osaka-city,
OSAKA, 530-8288 JAPAN

TEL No.: (+81)-(0) 6-6226-5256

FAX No.: (+81)-(0) 6-6226-5143

Authorized representative in the European Community:

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TEL No.: (+32)-(0) 16 381 211

FAX No.: (+32)-(0) 16 400 249

All brand names are trademarks or registered trademarks of TERUMO CORPORATION and their respective owners.

Table 1. Metacross RX Balloon Compliance Chart

Pressure (atm)	Balloon Diameter (mm)								
	3.00	4.00	5.00	6.00	7.00	8.00	9.00	10.00	12.00
8				5.97*	6.98*	8.06*	9.00*	9.91*	12.03*
10	3.02*	4.01*	5.04*	6.12	7.15	8.26	9.19	10.17	12.31
12	3.09	4.08	5.15	6.25	7.31	8.46	9.36**	10.44**	12.64**
14	3.16	4.15	5.26	6.40*	7.49*	8.66**	9.56	10.70	13.01
16	3.23	4.22	5.36	6.54	7.68	8.90			
18	3.30	4.30	5.48**						
20	3.38**	4.41**	5.61						
22	3.47	4.51	5.75						
24	3.57	4.64							

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© TERUMO CORPORATION September, 2021