

Navicross®

Support Catheters

MASTER YOUR PROCEDURES WITH PRECISE PERFORMANCE



NAVICROSS® 0.018"
Support Catheter

Now Available

 **TERUMO**
INTERVENTIONAL
SYSTEMS

PUSHING BOUNDARIES


Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.

We are relentlessly seeking new ways to help you achieve **better outcomes for more patients**.



NAVICROSS® 0.018" - Now Available

PRODUCT CODE	WIRE COMPATIBILITY	CATHETER LENGTH	TIP SHAPE
NC18650	0.018"	65 cm	Straight
NC18651	0.018"	65 cm	Angle
NC18900	0.018"	90 cm	Straight
NC18901	0.018"	90 cm	Angle
NC18130	0.018"	135 cm	Straight
NC18131	0.018"	135 cm	Angle
NC18150	0.018"	150 cm	Straight
NC18151	0.018"	150 cm	Angle

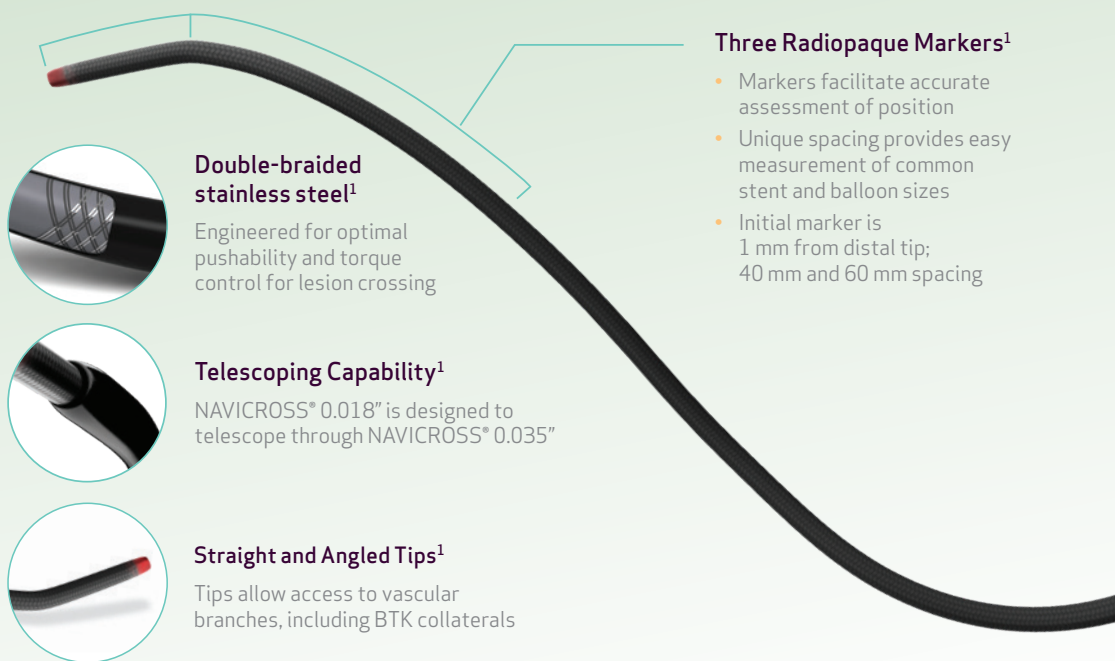

Double-braided stainless steel¹

Engineered for optimal pushability and torque control for lesion crossing

Telescoping Capability¹

NAVICROSS® 0.018" is designed to telescope through NAVICROSS® 0.035"

Straight and Angled Tips¹

Tips allow access to vascular branches, including BTK collaterals

Three Radiopaque Markers¹

- Markers facilitate accurate assessment of position
- Unique spacing provides easy measurement of common stent and balloon sizes
- Initial marker is 1 mm from distal tip; 40 mm and 60 mm spacing

FIND OUT MORE  Phone: 800.862.4143  terumo.com  Fax: 800.411.5870

RX ONLY. Refer to the product labels and package insert for complete warnings, precautions, potential complications, and instructions for use.

Reference:

1. Data on file.

©2019 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of Terumo. PM-02687

 **TERUMO**
INTERVENTIONAL
SYSTEMS