

PriorityOne[®] AC

Aspiration Catheter

**BETTER
DELIVERABILITY**
FOR EFFECTIVE
THROMBUS
REMOVAL



TERUMO
INTERVENTIONAL
SYSTEMS



PUSHING BOUNDARIES

Terumo Interventional Systems is **committed to your success** with innovative procedural solutions and ongoing support for your most challenging cases.

We are relentlessly seeking new ways to help you apply effective solutions and achieve **better outcomes for more patients.**

**TERUMO**
INTERVENTIONAL
SYSTEMS

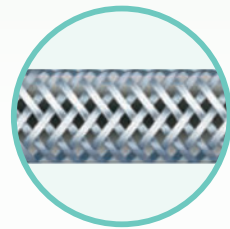
DELIVER EFFECTIVE ASPIRATION IN TORTUOUS ANATOMY

Designed to provide optimal deliverability for coronary and peripheral applications



KINK RESISTANCE for improved pushability

Fully braided stainless steel shaft provides variable stiffness with a stiffer proximal end for added pushability.

Guidewire-style stylet improves pushability and superior kink resistance—for optimal crossability in your most demanding cases.



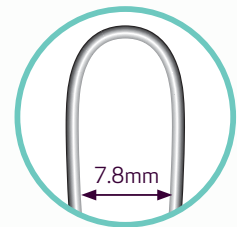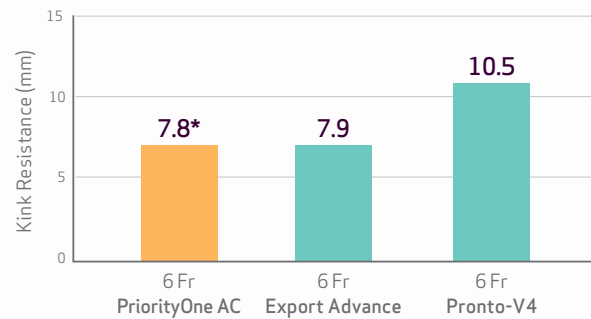
Proven Superiority in Distal Kink Resistance¹ Distance in millimeters when catheter kinked





*1.3x the kink resistance, compared to select competitors

SUPERIOR FLEXIBILITY for maximum trackability

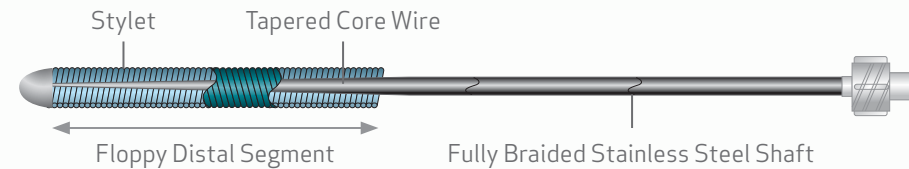Floppy distal segment and distal hydrophilic coating facilitates excellent trackability with superior flexibility.



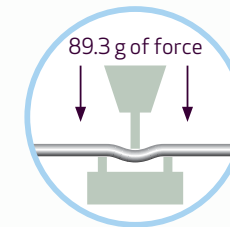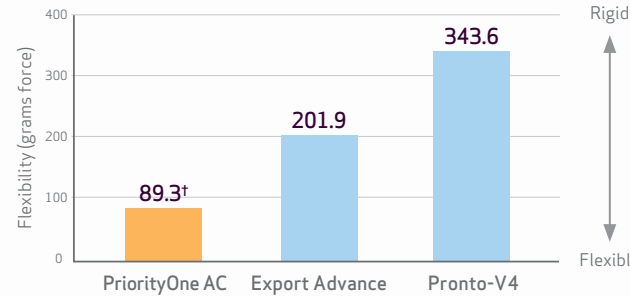Proven Superiority in Distal Flexibility¹ Force required to flex distal tip




[†]13x the flexibility, compared to select competitors

UNIQUE DESIGN for effective aspiration



Large extraction lumen reinforced by fully braided stainless steel shaft construction supports improved aspiration performance throughout the procedure.

Pre-loaded stylet and pre-connected extension line and syringe allow for easy setup.

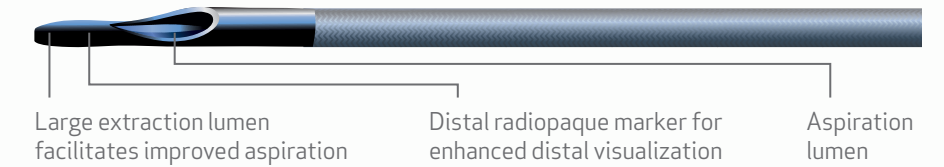Unique beveled tip design supports effective aspiration and supports deliverability.



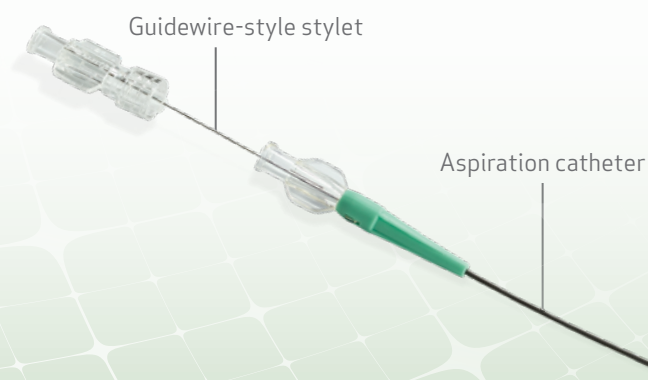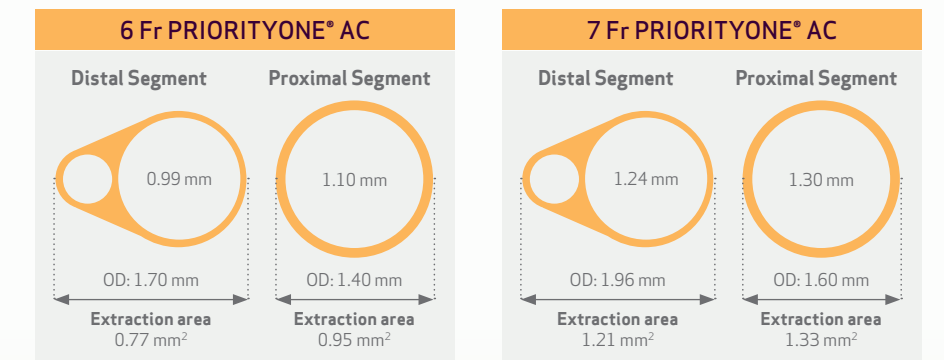Two sizes enable aspiration for a wide range of coronary and peripheral clinical scenarios.



PriorityOne® AC

Aspiration Catheter

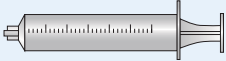

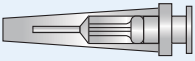
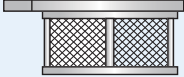
ORDERING INFORMATION

PRIORITYONE® AC	
PRODUCT CODE	GUIDE CATHETER COMPATIBILITY
30-6100	6 Fr (I.D. \geq 0.070" / 1.78 mm)
30-7100	7 Fr (I.D. \geq 0.080" / 2.03 mm)

PRODUCT SPECIFICATIONS

DISTAL TIP TUBE LENGTH	COATING	RAPID EXCHANGE SEGMENT	RADIOPAQUE MARKER BAND	POSITIONING MARKS	WIRE COMPATIBILITY	CATHETER LENGTH
9 mm	Hydrophilic	23 cm	1 mm	10 cm single white mark located 90 cm from distal tip	0.014" (0.36 mm)	140 cm

PRIORITYONE® AC COMPONENTS

30 ML LOCKING SYRINGES (2) 	EXTENSION LINE WITH 2-WAY STOPCOCK (1) 
GUIDEWIRE LUMEN FLUSHING TOOL (1) 	FILTER BASKET (1) 

FIND OUT MORE



US: 800.888.3786
CA: 833.883.7866



US: terumo.com
CA: terumocanada.ca

Reference:

1. Data on file, Terumo Clinical Supply Co., Ltd., Report No. 2025-TVC-009-00 (dated 11/10/2025).

RX ONLY. The advertisement is directed to physicians only, and not to consumers. Refer to product labels and packaging insert for complete warnings, precautions, potential complications, and instructions for use.

©2025 Terumo Medical Corporation. All rights reserved. All brand names are trademarks or registered trademarks of their respective owners. PM-05608 v2

**TERUMO**
INTERVENTIONAL
SYSTEMS